




Ascend[®] Full Face

replaceable AIR[°]gel cushion masks
featuring Advanced Cushion Technology

Rise to a higher level of comfort and convenience with the Ascend full face mask: our first-ever full face mask featuring a replaceable cushion. Our patented AIR[°]gel replaceable cushion is easily detached and reattached for easy cleaning and cushion replacement, saving both time and money.

Ascend to a Higher Level of Comfort and Convenience

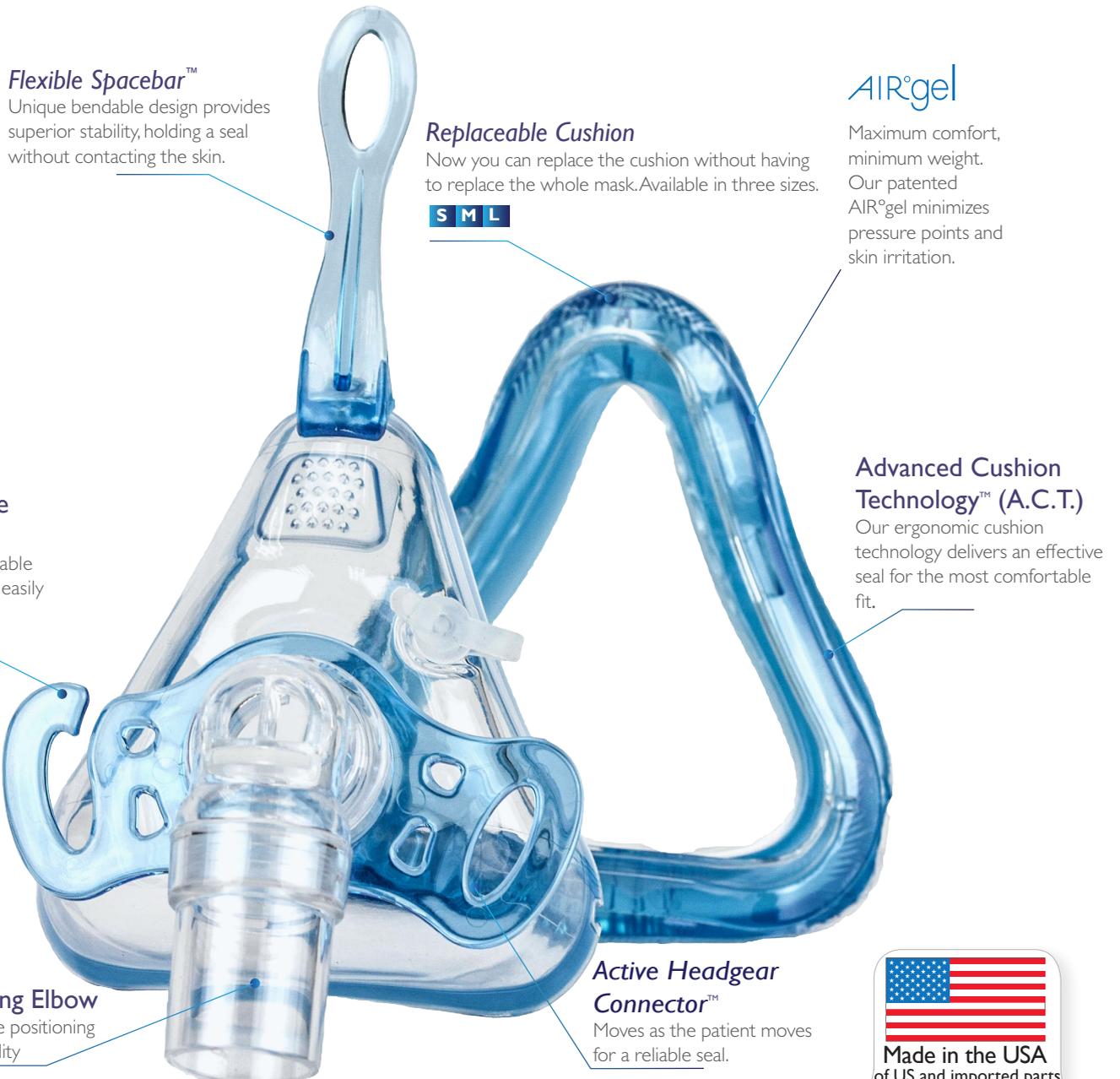

Flexible Spacebar™
Unique bendable design provides superior stability, holding a seal without contacting the skin.

Replaceable Cushion
Now you can replace the cushion without having to replace the whole mask. Available in three sizes.

S M L

AIR[°]gel
Maximum comfort, minimum weight. Our patented AIR[°]gel minimizes pressure points and skin irritation.

Quick-Release Headgear
Lightweight, breathable headgear slides off easily with one hand.

Advanced Cushion Technology™ (A.C.T.)
Our ergonomic cushion technology delivers an effective seal for the most comfortable fit.

360° Rotating Elbow
For optimal tube positioning and easier mobility

Active Headgear Connector™
Moves as the patient moves for a reliable seal.



What makes our Ascend Full Face Mask different?

- Our easy-to-use replaceable cushion allows you to extend the life of your mask.
- *Flexible Spacebar* is designed to provide superior stability without contacting the skin.
- *Active Headgear Connector* moves when the patient moves to help maintain an effective seal.
- Available in three sizes with simple and breathable quick-release EZFit headgear.

Advanced Cushion Technology (A.C.T.)

Our patented AIR°gel with an ergonomic design and luxurious feel reduces pressure while maintaining an effective seal.



Easy-to-Replace Cushion

Disassembly: Just squeeze the sides of the cushion in to release and snap in to replace.



Assembly: Place the top alignment pin in the notch and push the bottom ones in until the pins fit securely.



Ascend Full Face Masks with EZFit Disposable Headgear Specifications:

Product Number	Component	Leak Port	Entrainment Valve	Size	Mask Weight with cushion	Mask Dimensions with cushion
50825	Mask System & headgear	✓	✓	S, M, & L	NA	NA
50813	Mask & headgear	✓	✓	Small	94.5g	3.92"D x 3.86"W x 7.00"H
50814	Mask & headgear	✓	✓	Medium	103.4g	4.16"D x 4.20"W x 7.00"H
50815	Mask & headgear	✓	✓	Large	107.8g	4.24"D x 4.30"W x 7.00"H
50207	AIR°gel Cushion	✓	✓	Small	94.5g	3.92"D x 3.86"W x 7.00"H
50208	AIR°gel Cushion	✓	✓	Medium	103.4g	4.16"D x 4.20"W x 7.00"H
50209	AIR°gel Cushion	✓	✓	Large	107.8g	4.24"D x 4.30"W x 7.00"H
55051	EZFit Headgear	NA	NA	NA	NA	NA

CE 0123  Does not contain natural rubber latex

Ascend® Nasal
replaceable AIR°gel cushion masks

Also available as a Nasal Mask



Contact us today.

To place your order, call an authorized distributor.

Call 1.800.742.3646

Visit www.sleepnetmasks.com to learn more.

Sleepnet Corporation, 5 Merrill Industrial Drive,
Hampton, New Hampshire 03842 USA



Scan code to view Ascend Full Face Fitting Tips video.